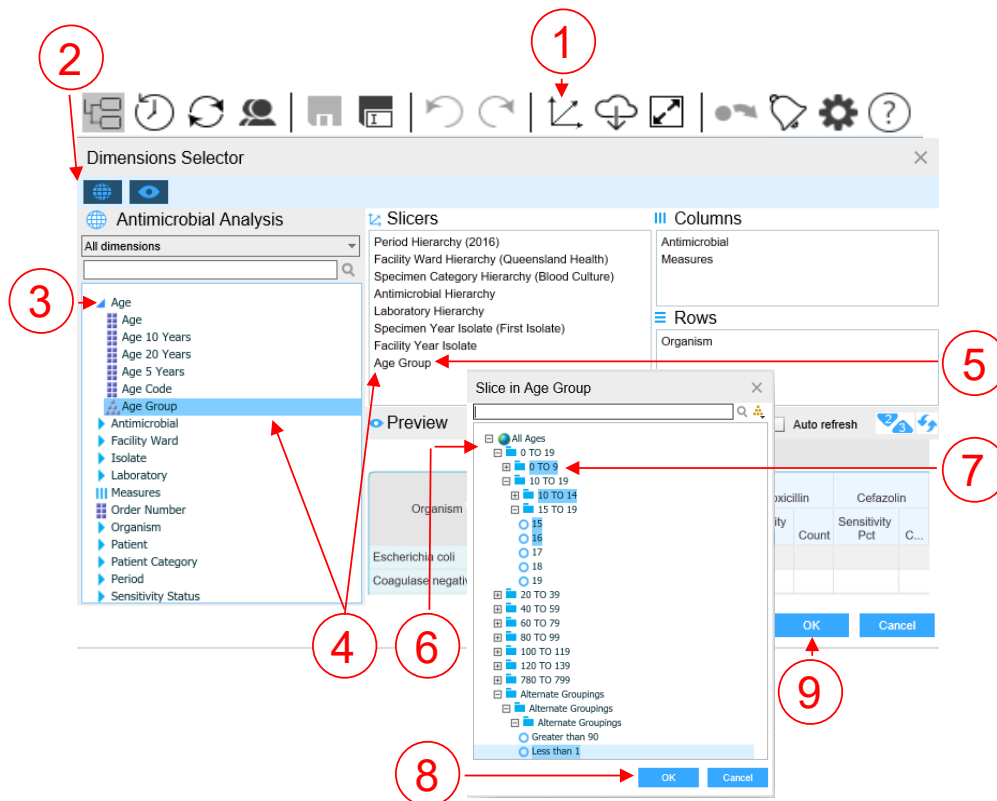









TRIM – D17-44650

OrgTRx Quick Reference Guide – Creating Paediatric Antibigrams

This document takes you through the steps to create a Paediatric Antibiogram. Open one of the Antibiogram views e.g. Antibiogram Blood.



- 1 Click the  *Dimension Selector*
- 2 Click  *All Dimensions* if necessary
- 3 Click *Age* dimension to expand
- 4 Click and drag *Age Group* to the *Slicers*

- 5 Double mouse click on Age Group to reveal the drop down menu with 20 year groups
- 6 Click on the  *Plus* sign to expand the groups to the required levels
- 7 Hold down the **Control key** and select the age range that meets your definition of 'paediatric'
For example, 0 to 9, and 10 to 14, and 15 and 16 and Less than 1
- 8 Click  *OK*
- 9 Click  *OK*

The result is a data set of isolates from children in the age range that has been selected.